


I need a hinge rail, do I need a double arm version with a drop down leg?

We offer a wide range of hinged rails for toilet and shower installations. There are single arm and double arm versions and a choice of rails with and without drop-down legs.

Here you must consider both the frequency of use and the likely size of the user(s).

Single Arm or Double Arm

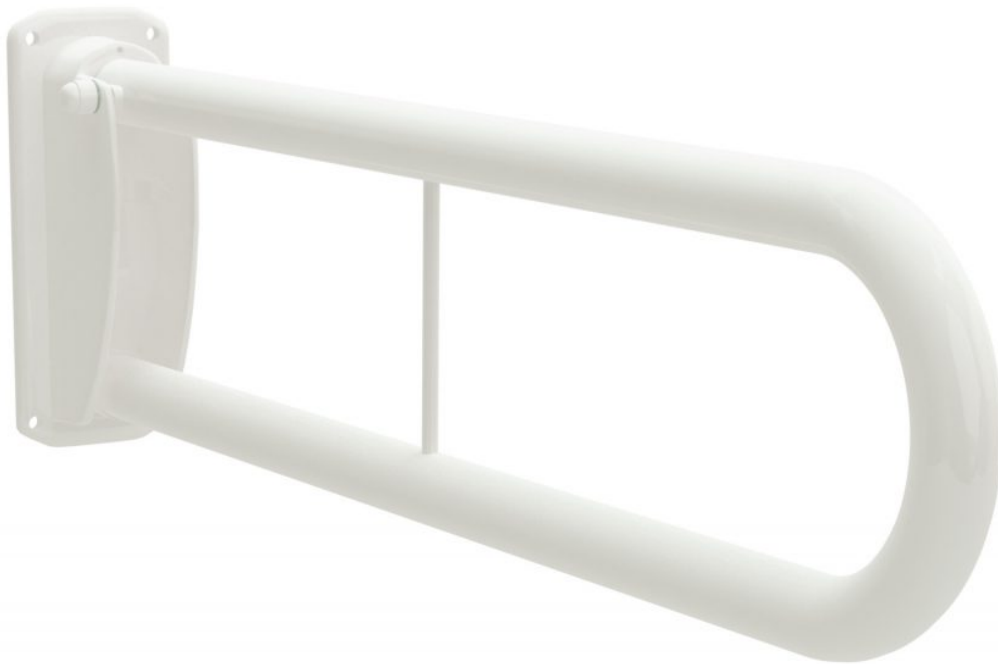

A single rail will provide stabilisation to your user up to a maximum loading of 130kg. A double rail will provide enhanced stability and will take a load bearing up to 170kg.

Drop Down Legs

A hinged rail with a drop down leg provides additional stability and security. The use of drop down legs does not increase the load bearing capability of the rails.

Friction or Slotting

Our friction hinged rails and shower seats improve the ease of use for those with more limited capabilities.

The friction hinge makes dropping and engaging the rail or seat easier with a smoother action. Slotted rails need to be lifted from a locked position to be lowered by the user.

